




**EverLift**<sup>®</sup>  
Submucosal Lifting Agent

**THE LIFT YOU'VE BEEN WAITING FOR**

# EverLift<sup>®</sup>

## Submucosal Lifting Agent

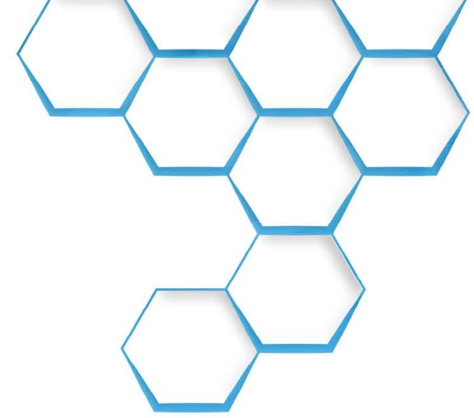

### INDICATION FOR USE<sup>2</sup>

EverLift Submucosal Lifting Agent is indicated for use in gastrointestinal endoscopic procedures for submucosal lift of polyps, adenomas, early-stage cancers, or other gastrointestinal lesions prior to excision with a snare or other appropriate endoscopic device.

**EverLift provides margin delineation with darker color.<sup>1</sup>**



### SUPERIOR CUSHION

#### CELLULOSE BASED & HYPERTONIC

EverLift is Cellulose based and able to stay in place longer because it does not get absorbed as fast as other lifting agents. It does not get absorbed as fast, because:



**Other Lifting Agents**

- EverLift is specially formulated with water retaining properties to prevent rapid absorption by the body.



**EverLift**

- EverLift is hypertonic, meaning our formulation has a slower rate of absorption because of osmotic resistance.

**EverLift holds a 2mm lift for 60 minutes.<sup>3</sup>**

# THE LIFT

you've been waiting for.

## CONVENIENTLY PACKAGED



### PRE-MIXED WITH METHYLENE BLUE<sup>2</sup>

Methylene blue is used to highlight the lesion perimeter and aid in: increased occurrence of en bloc resection rates, reduced occurrence of residual lesion and incomplete resections, improved identification of muscle injury to help avoid missed perforations.



### REQUIRES NO PRE-PREP

Convenient zero-prep pre-filled syringes attach directly to any endoscopic injection needle's luer-lock connection.

## COST EFFECTIVE

Use only what you need without throwing out what you paid for.

EverLift is reimbursable!  
Visit our reimbursement guide at [www.gi-supply.com](http://www.gi-supply.com).

First submucosal lifting agent designed to be cost effective and aid in complete resections.





**GI Supply has been the manufacturer of the leading endoscopic tattoos, Spot® and Spot® Ex, for 20 years and has been used in more than 5 million cases.**

GI Supply continues to expand its submucosal injection offering with innovative yet affordable devices to make your job easier. Founded by a gastroenterologist, we have your practice's best interests in mind.

ITEM NO.	DESCRIPTION	UNIT
<b>GIS-55</b>	EverLift Submucosal Lifting Agent 5 mL	10/Box
<b>GIS-59</b>	EverLift Submucosal Lifting Agent 10 mL	10/Box



**SUPERIOR CUSHION**



**CONVENIENTLY PACKAGED**



**COST EFFECTIVE**

- 1 - NAMSA Biological Safety Report. Evaluation of Carbon Black. July 2019
- 2 - G55-006 Rev 05 IFU
- 3 - Measurements taken by GI Supply. Data on file.

### Why choose EverLift?

"This is an effective lifting agent for endoscopic resection that is offered in convenient and cost-effective packaging."

#### Dr. Douglas Rex

Director of Endoscopy at  
Indiana University Hospital  
in Indianapolis



Scan to  
learn more!



**For More Information or Support**

Call: (800) 451-5797 Visit: [www.gi-supply.com](http://www.gi-supply.com)

©2021 GI Supply. All rights reserved. EverLift and Spot are registered trademarks of GI Supply Inc.

G-1587-04

